


# 014050 RRDE-3A/CS-3A ガスパージチューブ補修キット

弊社製品をお買い上げいただき有難うございます。  
本製品はRRDE-3A 回転リングディスク電極装置、CS-3A セルスタンド 用の配管チューブの補修キットです。

下記の弊社webサイトでチューブ交換の詳細な手順をご確認の上で作業を行ってください。

<https://www.bas.co.jp/2203.html#defaultTab12>


## 内容

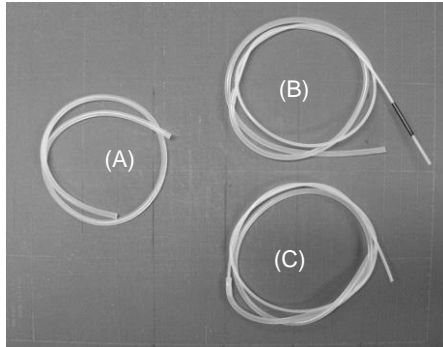

- (A) シリコンチューブ 1本
- (B) シリコンチューブ+テフロンチューブ (黒いカバー付き) 1本
- (C) シリコンチューブ+テフロンチューブ 1本

※取付作業を行うときは電源をOFFにし、電源ケーブルを取り外した状態で行ってください。

## RRDE-3A の配管チューブの交換

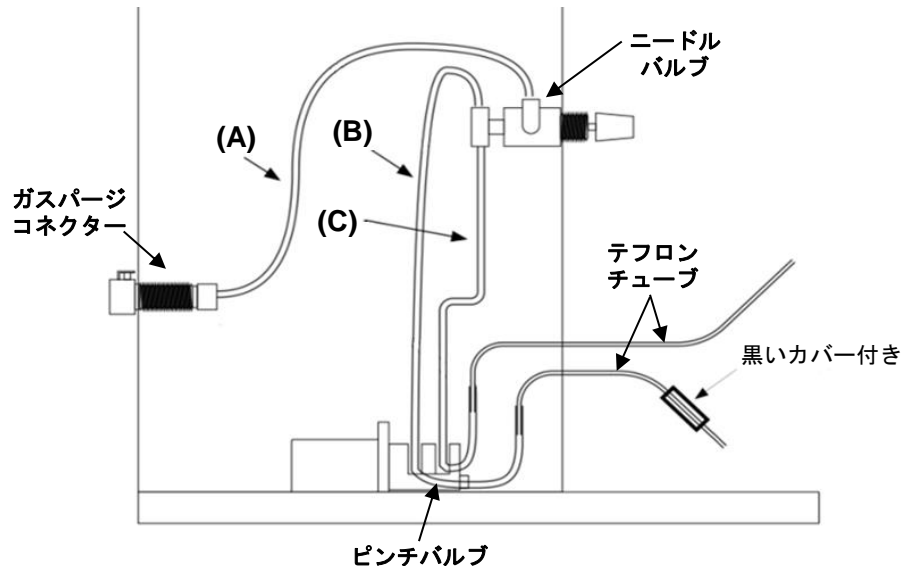


←背面

→前面

横から見た図

## CS-3A の配管チューブの交換

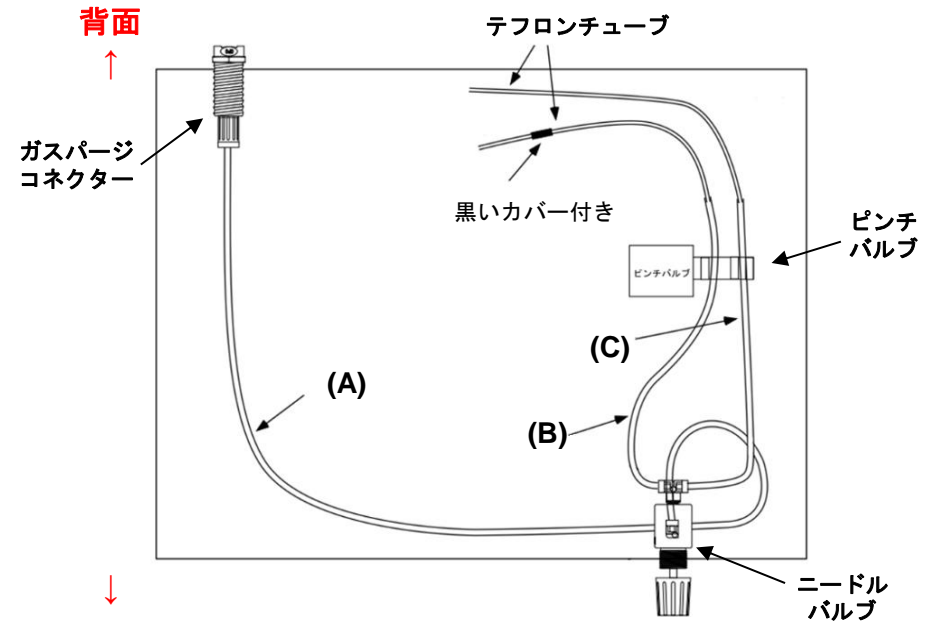

上から見た図

- (A) : ガスパージコネクタからニードルバルブまでのチューブ
- (B) : ニードルバルブから、ピンチバルブを通してパージ用チューブ (黒いカバー付き)
- (C) : ニードルバルブから、ピンチバルブを通してパージ用チューブ

### ・関連商品

- 012779 CS-3A セルスタンド Ver 1.1
- 013725 RRDE-3A 回転リングディスク電極装置 Ver.3.0
- 旧製品 RRDE-3A 回転リングディスク電極装置 Ver.2.0

ビー・エー・エス株式会社


<https://www.bas.co.jp>

E-mail: [sales@bas.co.jp](mailto:sales@bas.co.jp)

# 014050 RRDE-3A/CS-3A gas purge tube repair kit

Thank you for your purchase from BAS Inc.  
 This product is a tube repair kit for RRDE-3A Rotating Ring Disk Electrode Apparatus or CS-3A Cell Stand.  
 Please check the detailed instructions for tube replacement on our website below before proceeding.  
<https://www.als-japan.com/support-instrument-manual.html>


## Contents

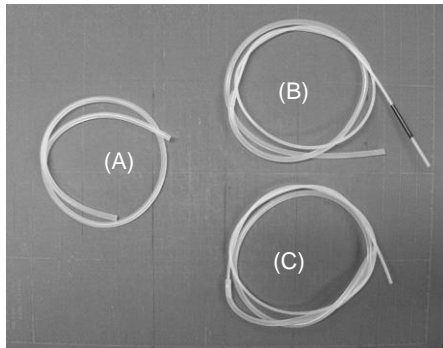


- (A) Silicone tube 1 pc
- (B) Silicone tube + Teflon tube (with black cover) 1 pc
- (C) Silicone tube + Teflon tube 1 pc

※Please turn off the power and unplug the power cable before performing replacement.

### RRDE-3A piping tube replacement

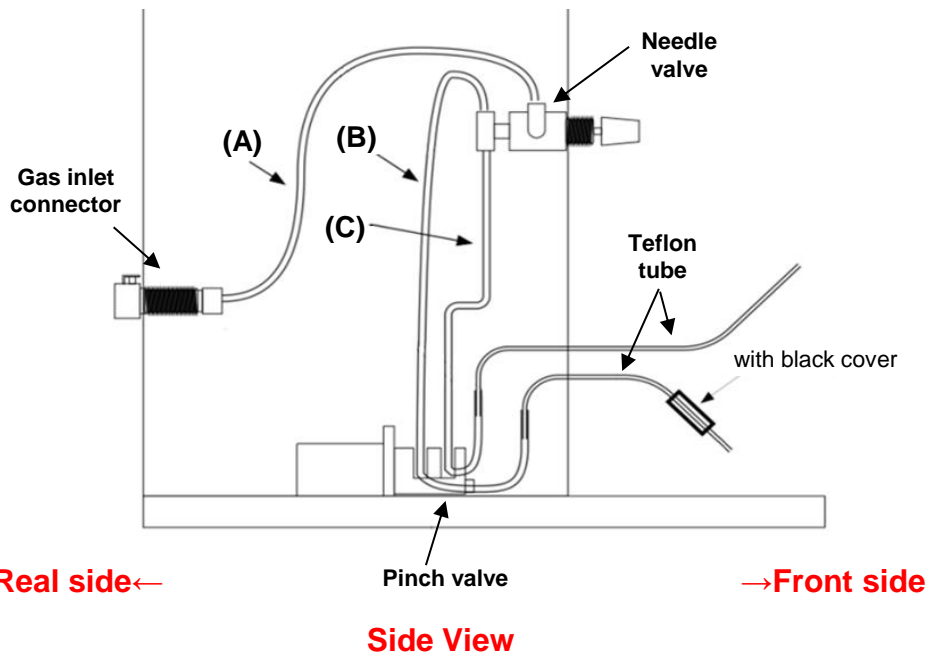

### CS-3A piping tube replacement

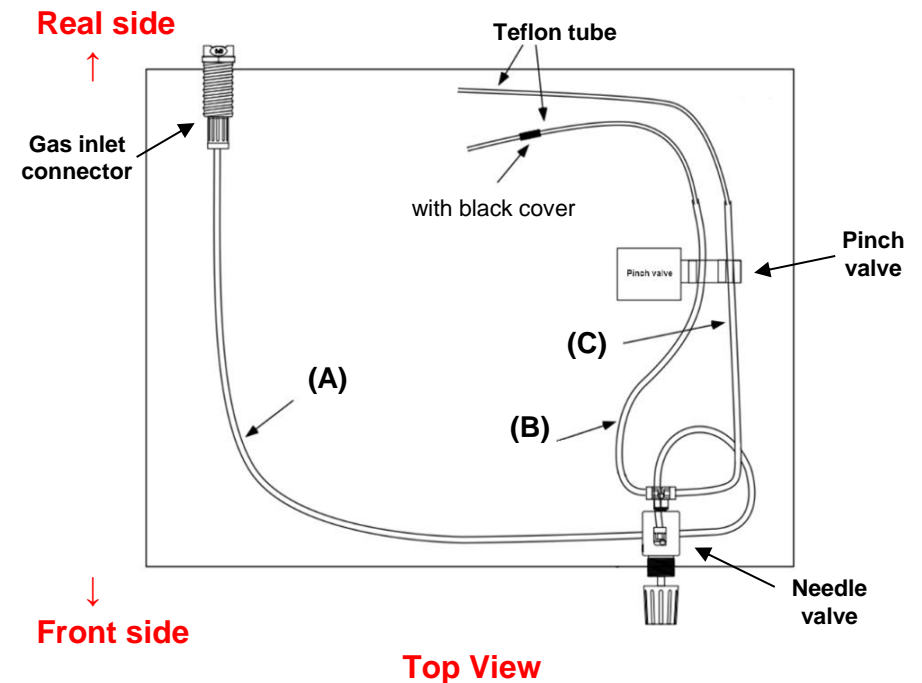

- (A): Tube from gas inlet connector to needle valve
- (B): Purge tubes from the needle valve through the pinch valve (with black cover)
- (C): Purge tubes from the needle valve through the pinch valve

- **Relative Product**
- 012779 CS-3A Cell Stand Ver.1.1
- 013725 RRDE-3A Rotating Ring Disk Electrode Apparatus Ver.3.0
- Previous product
- RRDE-3A Rotating Ring Disk Electrode Apparatus Ver.2.0